


Installation, Operation, and Maintenance Manual

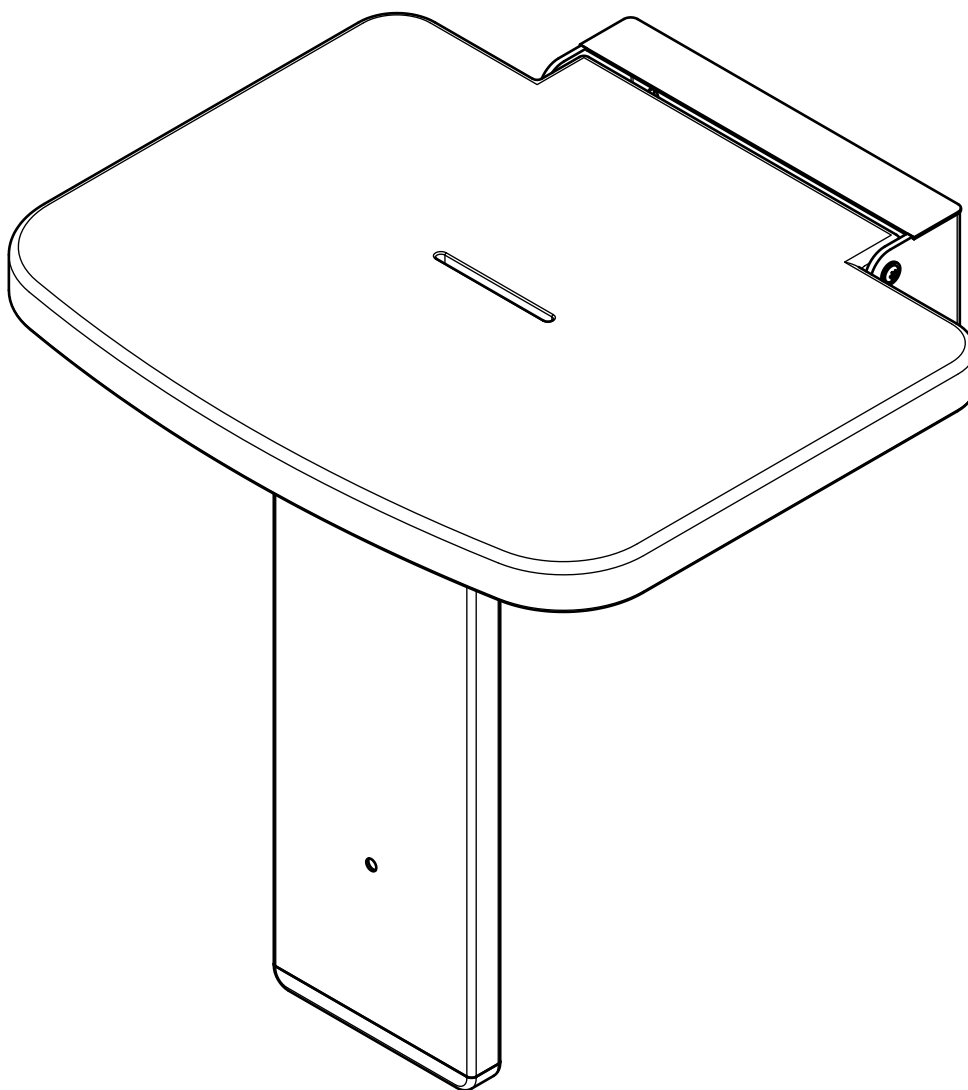

Product code: F22-FDS-SST

Description: FOLD DOWN SHOWER SEAT - PU PAD, STAINLESS STEEL

Disclaimer:

Fitzroy Of London pursues a continuous improvement strategy in the design and performance of its products and reserves the right to change specifications without prior notice - terms and conditions apply. The information contained in this document is correct at date of issue. Dimensions are provided as a guide only - some variation may occur due to manufacturing tolerances.

Fitzroy Of London accepts no liability whatsoever for personal injury or damage to product(s) or property that has been caused by incorrect assembly, installation, modification or use of this product. Failure to install/use this product in line with the instructions in this manual will invalidate any warranty and may affect product safety.

Safety warnings:



Always read the installation, operation, and maintenance manual. Installation should be carried out by a suitably qualified tradesperson using the correct wall fixings for the wall substrate. Familiarise yourself with these warnings before beginning any work:

- Installers: leave a copy of this manual with the end user (printed, digital, or access to QR link)
- Pre-drilling fixing holes is not recommended before physically receiving this product.
- Never alter the fixing holes and use all fixing holes provided.
- A thorough check must be carried out upon receipt and any defects reported immediately.
- Opening the packaging without due care and attention could result in damage to the product or injury.
- This product must be stored in a dry place, the correct way up and not stacked under other equipment.
- Check there is no risk of contacting/damaging electrical wires or water/gas pipes before attempting to drill into the wall.
- Ensure that the surface which the product is to be fitted is flat and free from defects (wall and floor).
- Ensure any shower trays/raised shower platforms are suitably supported or firm enough to withstand the pressure from the feet of the product.
- Non-corrosive fixings must be used to fix this product to the wall - **supplied wall plugs are designed for solid block walls only, NOT dry walls - Fitzroy of London recommends product is fitted to solid block walls.**
- The product(s) should be installed by a suitably competent person or qualified tradesperson.
- The installer must thoroughly check the integrity of the assembled and installed product.
- Improper installation/use may cause accidents and injuries.
- Modification to the product is not recommended and may result in accident / injury and invalidates your warranty.
- This manual should be retained by the end user/property owner for the life of the product.
- The product should not be subjected to impact.
- Do not stand on the installed product.
- Do not use the product if it appears damaged. Check before each use.
- Do not smoke near, or when using, the product. Do not expose to naked flames.
- If in doubt, carry out your own risk assessment and/or quality check for each installation - not recommended for high-risk settings.
- Recommend removal from setting if signs of misuse are present.
- Stand clear when lowering seat. Raise to stowed (vertical) position after each use.
- **Maximum load = 200 kg**

Product cleaning:



CAUTION - Improper cleaning can lead to surface damage and injury.

- For the avoidance of bacterial/mould build-up, regular deep cleaning is recommended, paying attention to all surfaces, crevices, joints and recessed areas on the product.
- After each use, rinse the seat with clean water to remove soap and water traces - wipe dry with a soft cloth.
- Follow the cleaning product manufacturer's advice about contact time, concentration, safety etc.
- Try on an inconspicuous area first.
- Rinse with fresh water immediately after cleaning.
- Do not use any cleaning products with an unknown composition.
- Do not mix chemicals.
- Note the material specific information on suitable/unsuitable cleaning products.
- Thorough drying with a soft dry cloth is recommended.

Unsuitable cleaning products:





- Abrasive products (e.g., scoured/scouring agents, brushes, steel wool, dirt easers, microfibers).
- Acidic, alkaline, chlorine and solvent based cleaning products.
- Steam cleaners.
- **DO NOT** use cleaning components containing iron (e.g., scour sponge, steel wool).

NOTE: abrasives and harsh chemicals are especially damaging to the PU (polyurethane) seat pad.

Disposal:

Dispose of this product and packaging in line with local authority regulations - deconstruct if necessary.

Product dimensions: (mm)

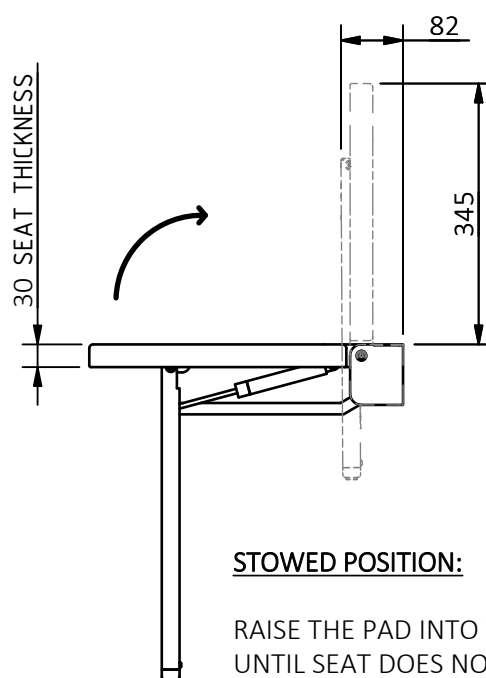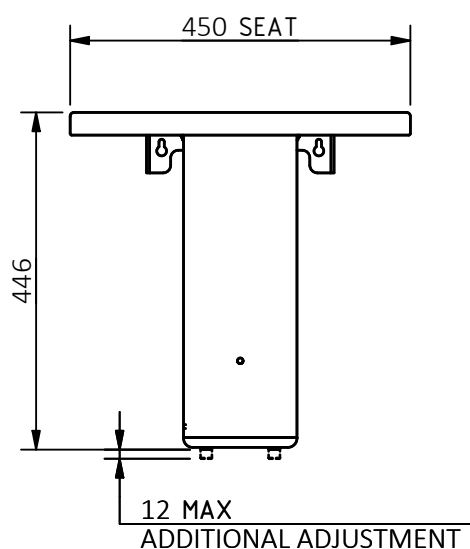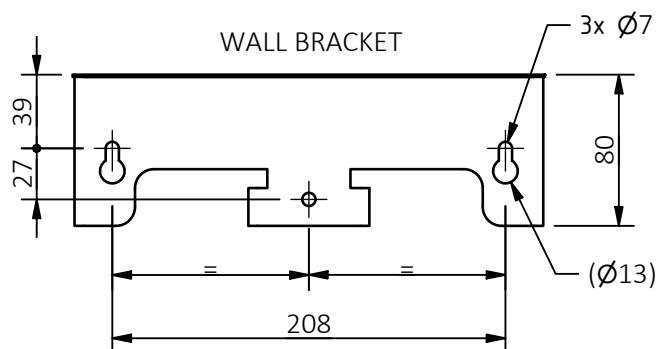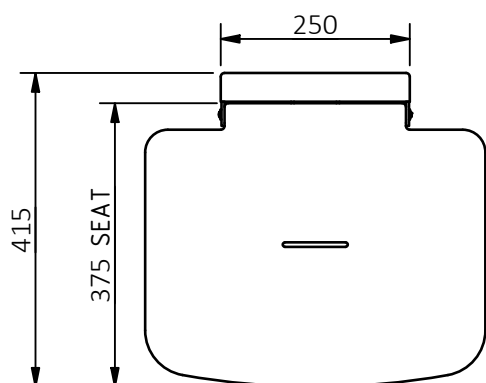

STOWED POSITION:


RAISE THE PAD INTO VERTICAL POSITION UNTIL SEAT DOES NOT DESCEND (VERTICAL OR JUST BEYOND VERTICAL).

GAS STRUT PROVIDES SLOW-DESCEND MECHANISM FOR PAD AND LEG.


STAND CLEAR WHEN LOWERING.

LEAVE IN STOWED POSITION AFTER EACH USE.

Tools required:

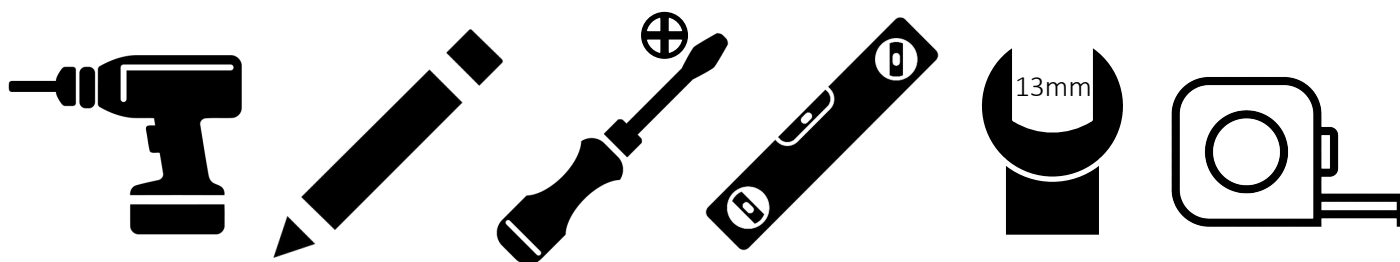

What's in the box:

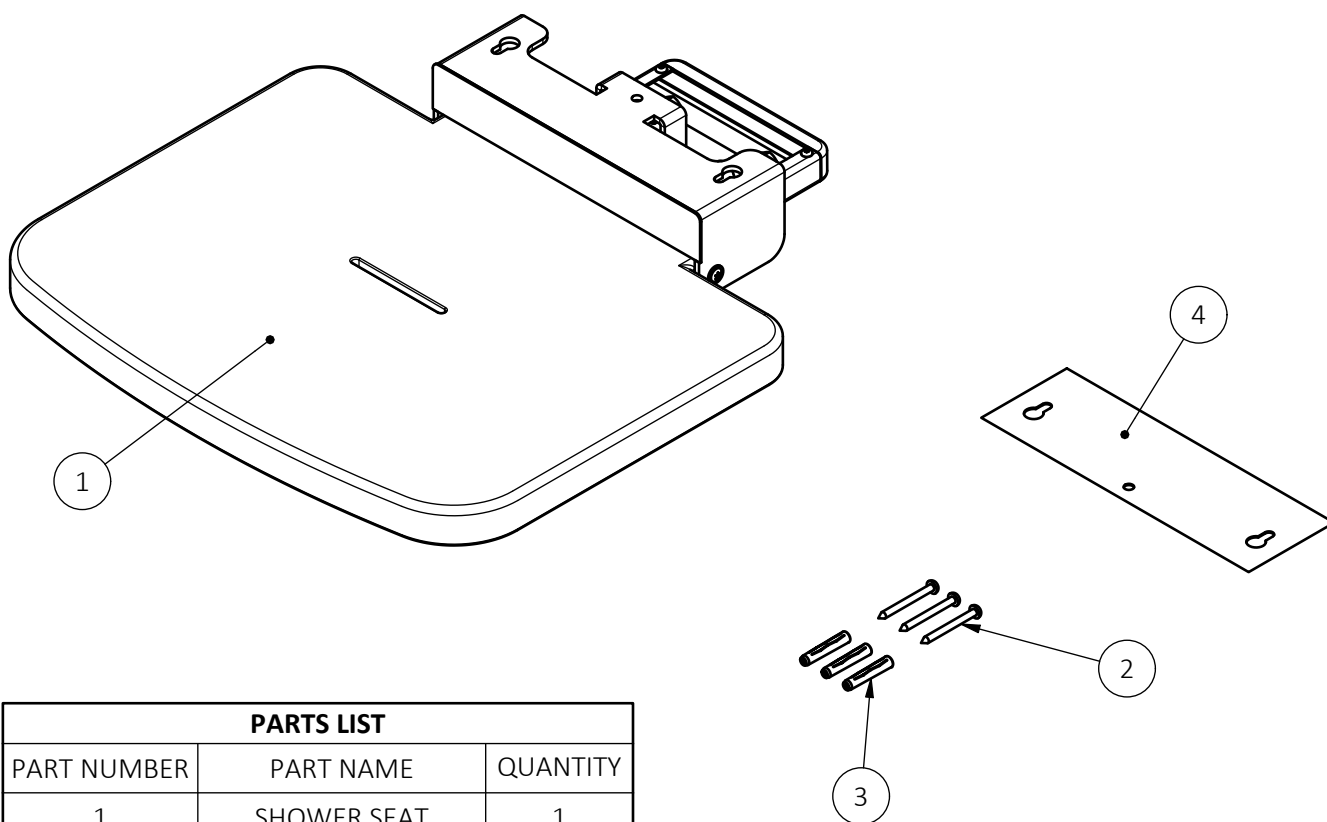

PARTS LIST		
PART NUMBER	PART NAME	QUANTITY
1	SHOWER SEAT	1
2	SCREW FIXINGS	3
3	PLUGS (FOR BLOCK WORK ONLY)	3
4	FIXING TEMPLATE	1

Pre-installation:




CAUTION - Ensure safety warnings section of this document has been read.

Installation:

1.

Remove all parts from the packaging provided.

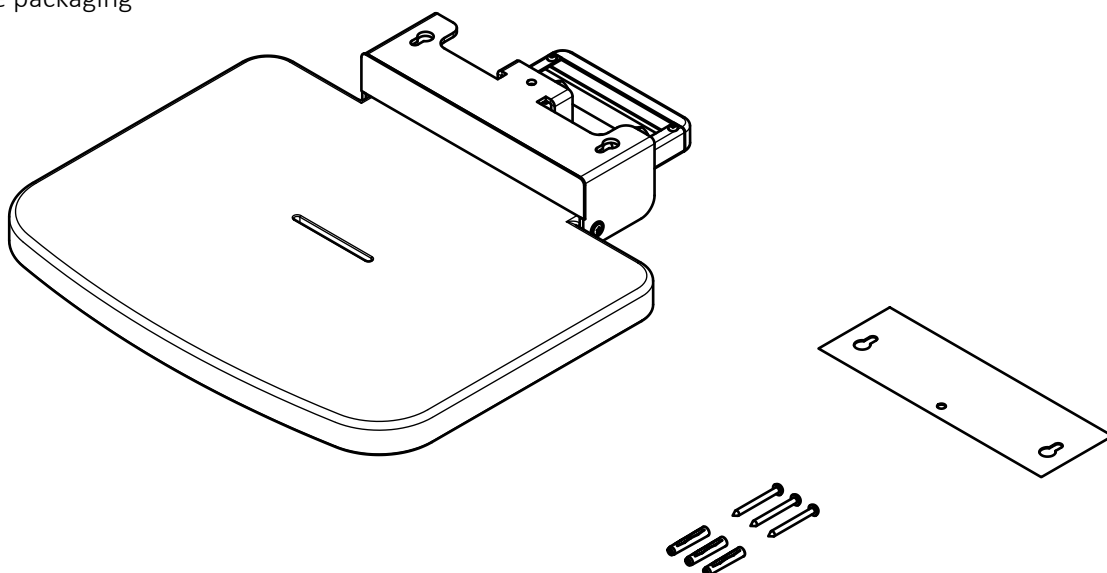


2.

Mark the required height at which the shower seat needs to be installed using the top edge of the wall bracket.

Ensure that the seat is level.



This step may be easier with 2 people: 1 to hold, 1 to mark the positions and check if level.

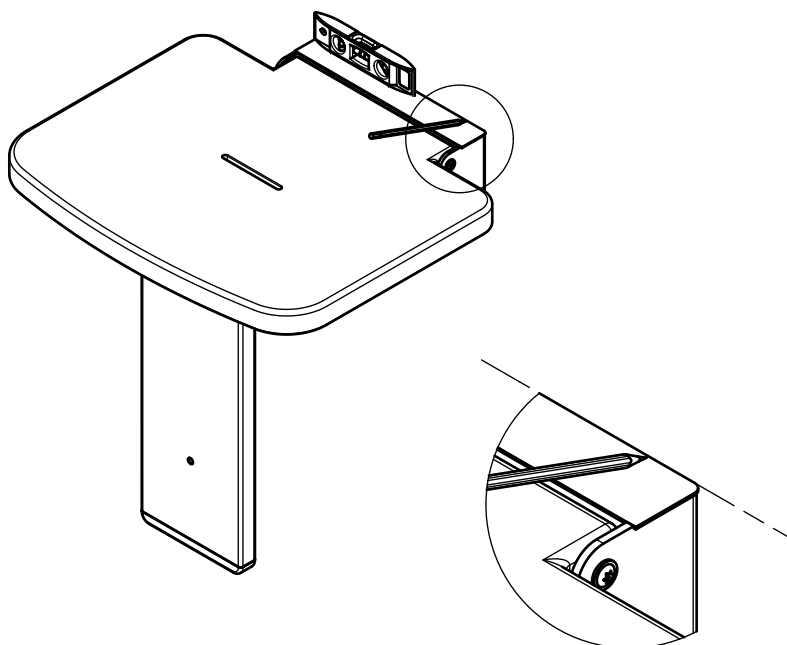


Installation:

3.

Peel off the back of the template.

Stick the installation template against the wall, ensuring it is level.

Make sure the top edge of the template is on the drawn line.

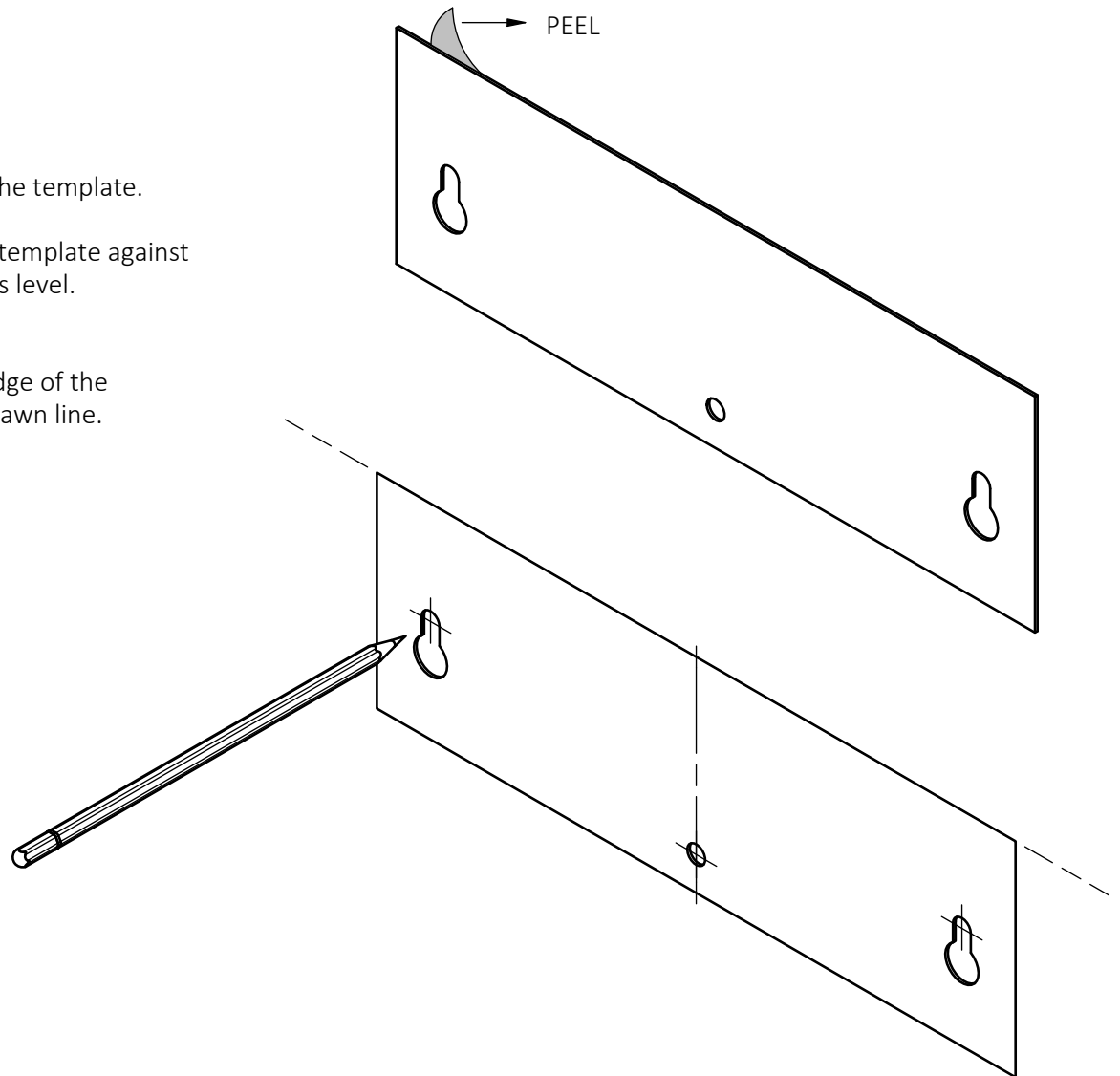


4.

Referencing the drawing on Page 4, and using the template adhered to the surface, drill the required fixing hole locations to the correct diameter and depth for the supplied ($\varnothing 8\text{mm}$) plugs.



Drill hole at the top of the keyholes as shown.

Remove the template sticker once complete.

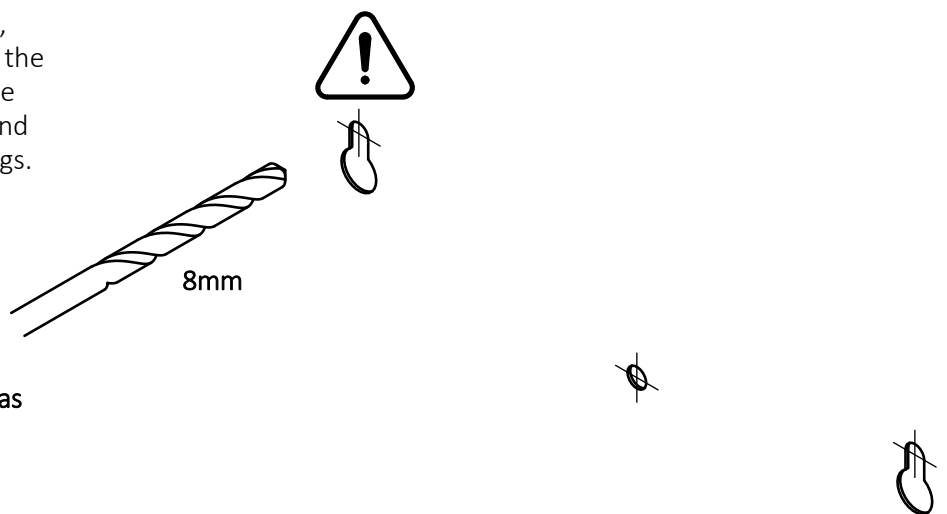

Installation:


5.

Insert plugs into ALL fixing holes.



ENSURE PLUG FULLY INSERTED INTO BLOCK WORK (PUSHED BEYOND TILE LAYER).

SUPPLIED PLUGS FOR DENSE BLOCK-WORK ONLY! (SEE SAFETY WARNINGS)

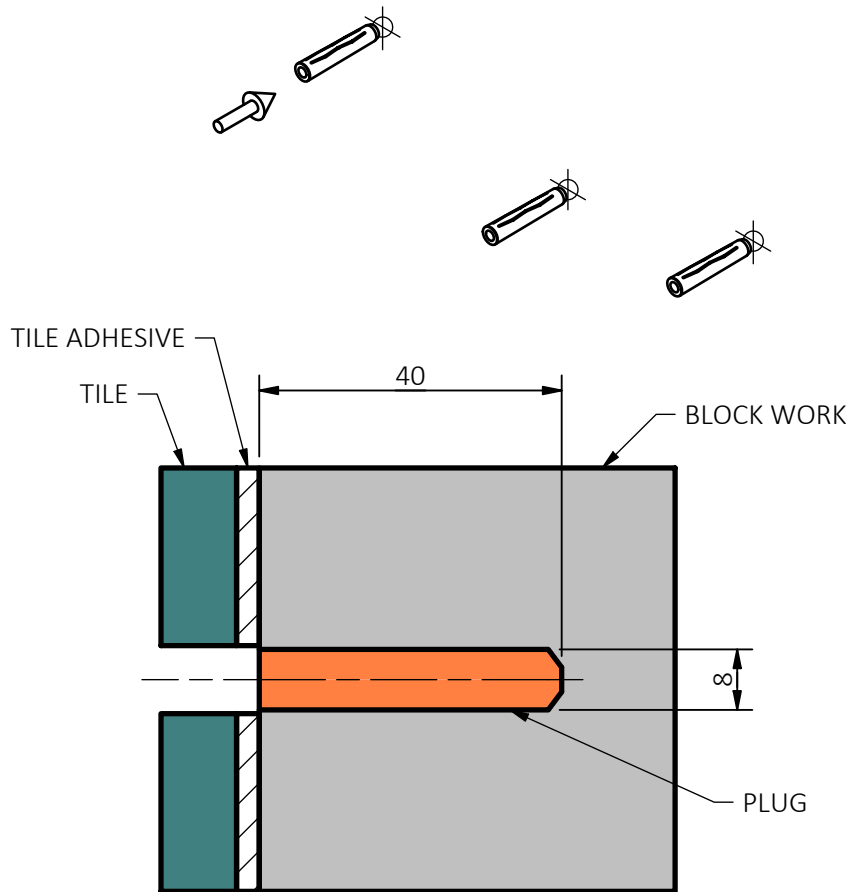

6.

Insert supplied screws into the **LEFT and RIGHT** fixing holes **only**.


DO NOT insert fixings completely - leave proud by at least 5-10mm.

Stainless fixings **MUST** be used.

Ensure all fixing heads are suitable to fit within the frame



Seal around ALL wall hole boundaries with suitable silicone sealant - highly recommended to reduce chance of water ingress!

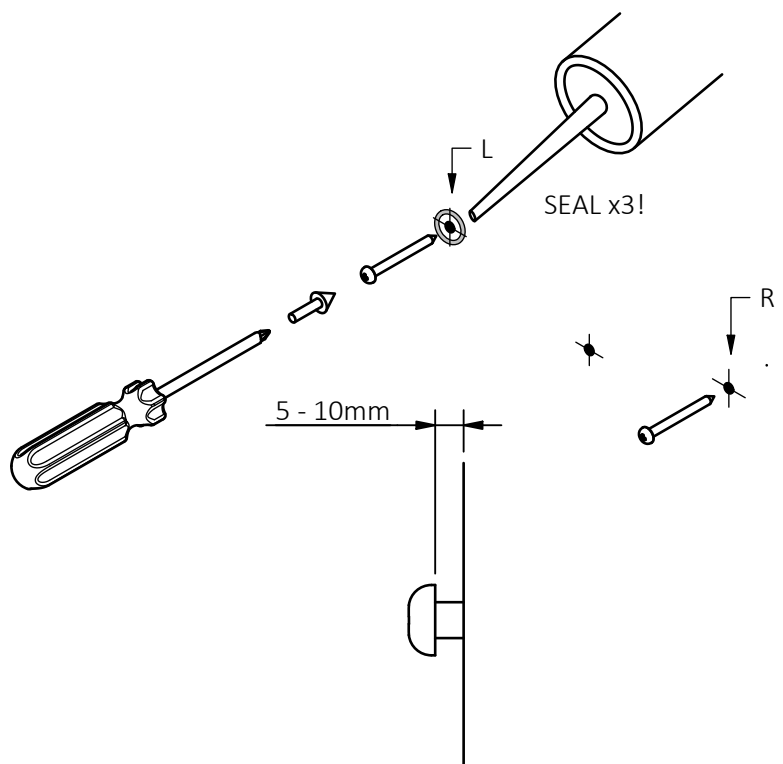


Installation:


7.

Locate the wall bracket over the LEFT HAND and RIGHT HAND screw fixing heads using the keyholes.

ENGAGE the wall bracket downward and into position until resting on the screws.

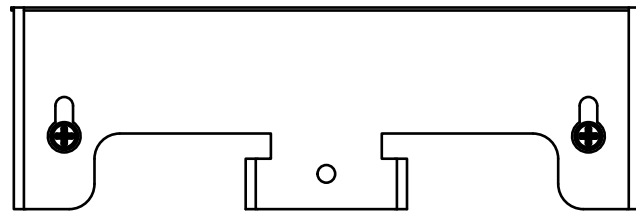




ENGAGE

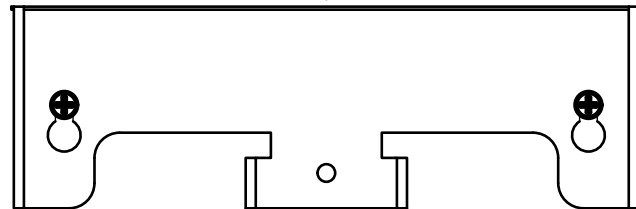

8.

Tighten the LEFT HAND and RIGHT HAND screws fully - (pressing against the wall bracket).

Add the final central screw to the centre hole - Tighten fully.

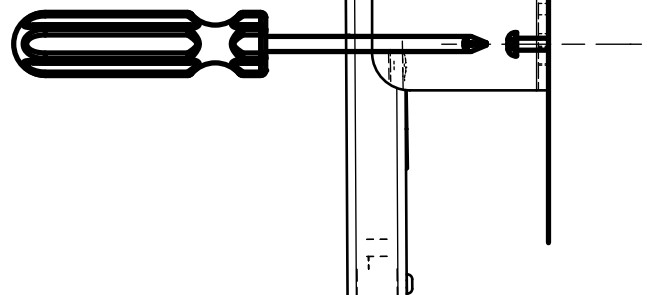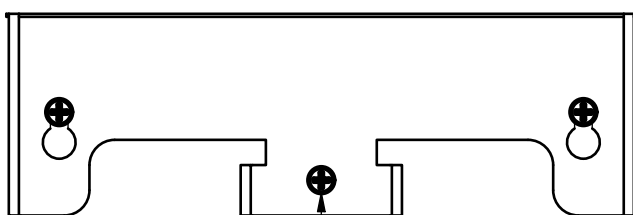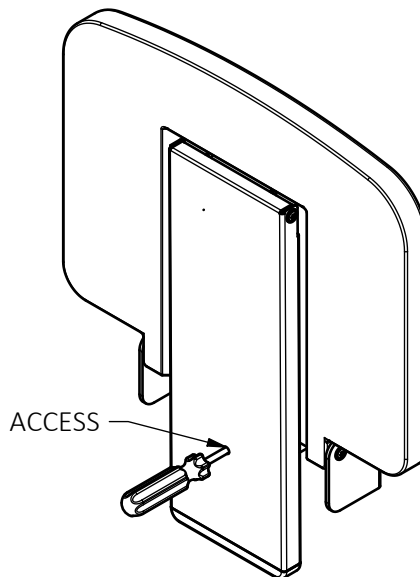


Installation:

9.

Feet Fine Adjustment

If the feet of the shower seat do not sit firmly on the ground - the feet can be adjusted manually.

1. Loosen the locking nut at the top of the feet.
2. Turn the rubber feet to adjust the level as necessary.
3. Re-tighten the locking nut at the top, securing the foot position at the required length.

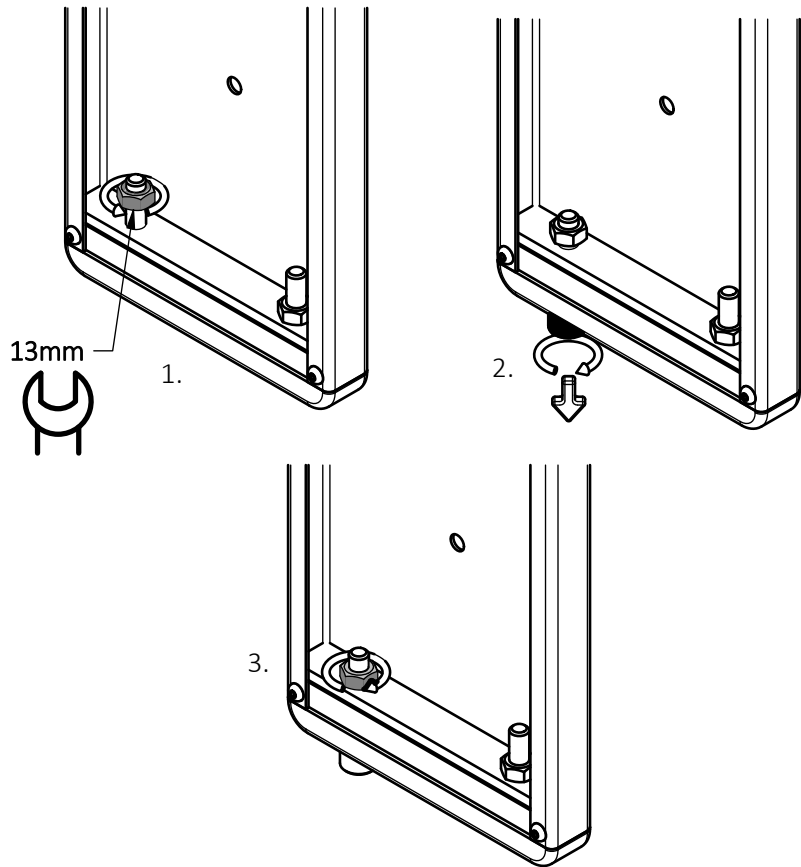


10.



Seal top edge and sides of wall bracket with suitable silicone sealant.

See also, STEP 6.

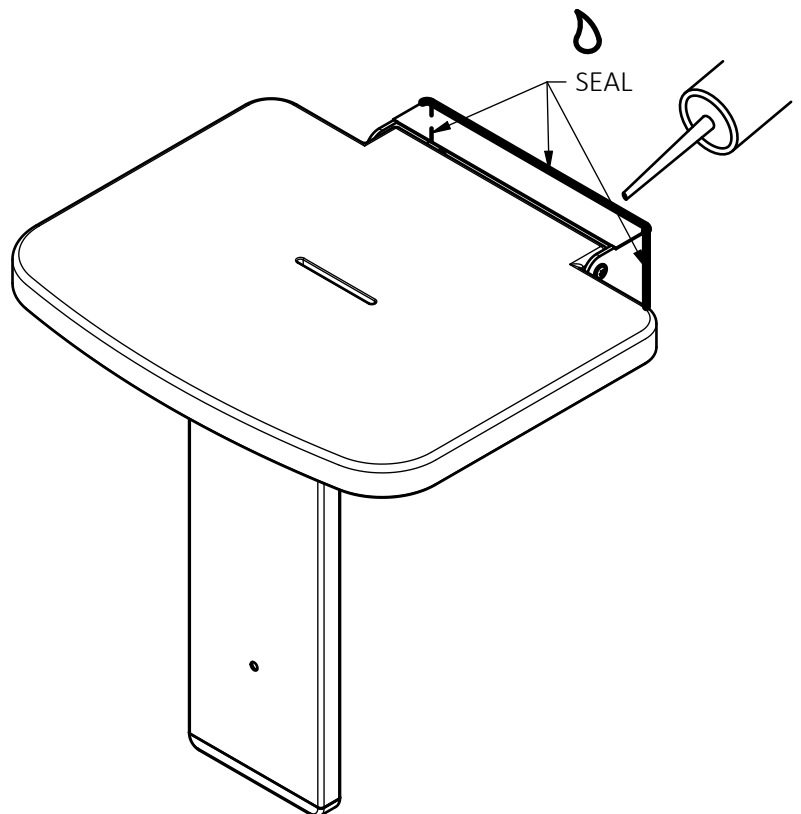

Installation:


11.

Perform a visual and physical check to ensure the positioning of the seat is correct and that it is safe to use.

Check up/down operation and level leg position, plus security on the wall.



Check shower tray/floor area does not flex excessively under load - additional support may be required under the tray area - see also Safety warnings section of this manual.

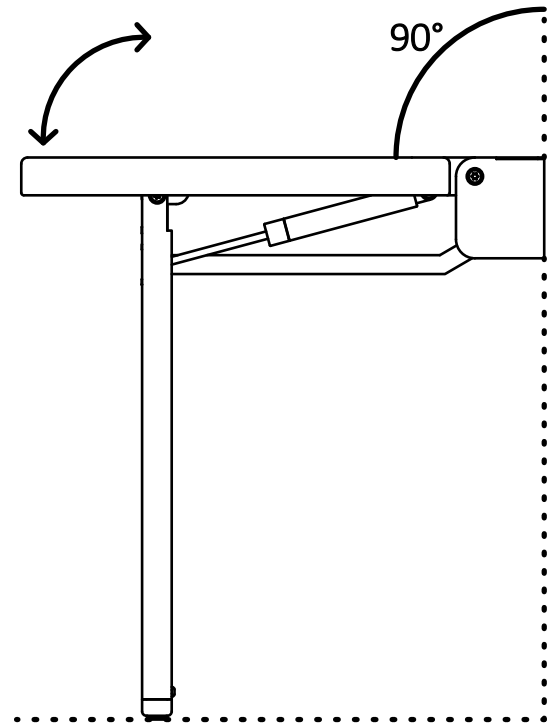

GENERAL OPERATION ADVICE:

RAISE/LOWER THE SEAT BY HAND GENTLY, USING THE FRONT-MOST EDGE OF THE SEAT PAD.
DO NOT FORCE THE MECHANISM UP/DOWN.

ALWAYS LEAVE IN THE STOWED POSITION AFTER EACH USE.

END OF DOCUMENT